

Reliance



Orthodontic Products



Reliance Orthodontic Products, Inc.

1540 West Thorndale Ave., PO Box 678, Itasca, IL 60143
www.relianceorthodontics.com

Dear Doctor ,

Welcome to the Reliance family of orthodontic adhesive products!

As you will notice, we specialize in adhesive for all bonding needs, with all types of appliances, to any surface. Since our inception in 1981, we have achieved the reputation of providing consistent, quality adhesives, in a variety of delivery systems, at a cost far below competition. Before we sell a product commercially, it is clinically field tested. This thorough evaluation insures that each product performs to satisfaction for the full treatment time and is chemically safe for you, your staff and your patients.

As an added service to our adhesive customers, any bonding questions you have can be answered by calling me toll free at 1-800-323-4348 or email me at pgange@relianceorthodontics.com. We sincerely appreciate your patronage and look forward to continuing to serve your requirements in the months and years to come.

Paul Gange
President



Reliance Orthodontic Products, Inc.

Corporate Headquarters
1540 West Thorndale Ave.
PO Box 678
Itasca, IL 60143

TERMS AND CONDITIONS

ORDERING

To order Reliance Orthodontic Products, simply contact our home office Toll Free (in U.S. and Canada) at 1-800-323-4348 or 1-630-773-3580. You can also order through our secure website at www.relianceorthodontics.com. There is a minimum of \$25.00 on all orders.

PAYMENT

Options: Net 30 days, upon approval. Past due balances may be charged 1.5% per month finance charge at the sole discretion of Reliance Orthodontics. We accept Mastercard, Visa and Discover for any purchase.

SHIPPING

All prices are F.O.B. Itasca, and are shipped via FedEx. All shipping is billed with the order. Alternate shipping services are available. Please call for more information.

EXPORT


Export orders must be paid for in advance by proforma invoice.

TERMS

Net 30 days. Illinois residents please add 1.75% sales tax. An invoice is enclosed with each shipment and payment should be made from this invoice. Any outstanding balance will be subject to a finance charge of 1.5% per month. All payments are due in U.S. funds

EXCHANGES or REFUNDS

Your satisfaction is Guaranteed! Reliance Orthodontic Products, Inc. guarantees that its products are of the highest quality and workmanship. If at any time prior to its expiration date your adhesive system fails to perform in a manner not consistent with the product, we will replace the remaining product with new material free of charge or issue a credit for the remaining product. No refunds, credits or replacement will be issued after 60 days unless the product is deemed defective by Reliance Orthodontic Products, Inc. Products that are being returned to Reliance Orthodontic Products, Inc. must have a prior Return Authorization Number (RMA.) All products being returned for credit or replacement must be sent postage paid with a copy of your invoice and your RMA number clearly visible on the package, or the package will be refused by Reliance Orthodontic Products, Inc.

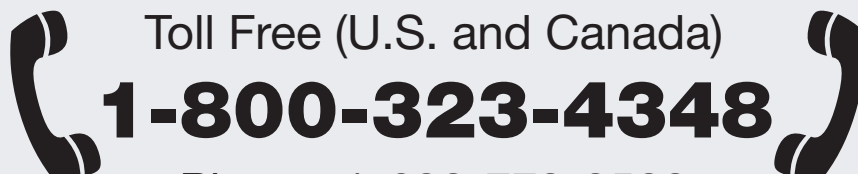

Toll Free (U.S. and Canada)

1-800-323-4348

Phone: 1-630-773-3580

Fax: 1-630-250-7704



Visit us online at

"www.relianceorthodontics.com"

Please email any questions or
comments directly to
pgange@relianceorthodontics.com

or

paulgangejr@relianceorthodontics.com



An ISO 13485:2003 Certified U.S.A.
and F.D.A. Registered Company

Reliance Products are proudly manufactured
and packaged in the U.S.A.

Please Note:

*We are constantly adding new products to keep abreast of new technologies.
If you do not see the item you need, or for details on future products,
please contact us at 1-800-323-4348.*

Chemical Cure Adhesives	Rely•a•Bond®	5
	Phase II®	6
Light Cure Adhesives	Light Bond™	7
	Pad Lock®	8
	Assure®	9
	Quick Cure™	9
Light Cure Accessories	UV Disclosing Light	10
	Uvex Skyper® Glasses	10
	Uvex Ultra Spec® 2000 Glasses	10
	Power Slot™	10
Enamel Protectant Sealants	PRO SEAL®	11
	L.E.D. PRO SEAL®	11
	Maximum Cure® Sealant.....	11
	Maximum Cure® Filled Sealant.....	11
Light Cure Retainer Adhesives	L.C.R.™	12
	FlowTain™	12
	FlowTain™ L.V.....	12
Specialty Adhesives	Excel®	13
Dual Cure Adhesives	Phase II® Dual Cure.....	14
Band/Crown Cements	Band-Lok®	15
	Crown Lok™	15
	Self Mixing Dual Syringes (Band-Lok® and Crown Lok™).....	15
	Ultra Band-Lok®	16
	Ultra Band-Lok® Luer Loc	16
	Precedent®	16
Ortho Veneer System	O.V.S.™	17
Indirect Bonding	Therma-Cure™	18
	Custom I.Q.™	18
Adhesion Boosters	Enhance™ Adhesion Booster	18
	Enhance™ L.C.....	18
	Assure® Plus	19
	Assure®	19
Enamel Conditioners	S.E.P.™ Self Etching Primer	20
	Etchant (Liquid and Gel)	20
	Blue Gel Etchant.....	20
Surface Conditioners	Porc-Etch™ (Porcelain Etchant).....	20
	Porcelain Conditioner	20
	Plastic Conditioner	20
	Metal Primer.....	20
	Rescue™ All Surface Bonding Kit.....	20
Prophylaxis Materials	1st & Final®	21
	Prep Step™ Propphy Angles.....	21
	Access™ Cheek Retractor.....	21
Finishing Systems	Renew® Finishing System (Burs and Points).....	22
	Restore™ Polishing Paste.....	22
Clinical Training Aids	Bonding Reference Guide	23
	Bonding Technique Cards.....	23
Lingual Retention and Kobayashi Ties	RETAINIUM™	24
	Stainless Steel Ortho FlexTech®	24
	Bond•a•Braid®	24
	Ortho FlexTech®	24
	Kobayashi Tie Hooks	24
	Pre-twisted Kobayashi Tie Hooks	24
Chairside Accessories	Brushes	25
	Sponge Pellets.....	25
	Spatulas.....	25
	Mixing Wells	25
	Mixing Pads	25
	Mark II CR Syringe.....	25
	EtchMaster® System	25
	EtchMaster® Regular and Lingual Tips	25
Aligner Adhesives	Bond Aligner™	26
Aligner Pontic Paint	Perfect•A•Smile™	26
A² Aligner Attachments	A ² ™ Aligner Attachment Kit.....	27
Oral Surgery Extrusion System	C&C™ Kit (Composite & Chain).....	27

Rely•a•Bond®

No Mix Adhesive



Contents of Kit

- 14 g Adhesive Paste (Jar, Syringe Tips or Push Syringes) •
- 16 g Primer, with or without fluoride •
- 9 g Etching Agent (Liquid or Gel) •
- Mixing Pad • Spatulas • Mixing Well •
- Brushes (Sponge pellets may be substituted for brushes) •

Product Features

- **Contains Fluoride** - Rely•a•Bond® Primer is formulated with a monomer that contains fluoride and will not reduce bond strength. Non-fluoride formula also available.
- **Unsurpassed Bond Strength** - Rely•a•Bond® is formulated with a highly filled paste that combines with a unique primer to produce a bond strength to etched enamel that virtually eliminates bonding failures.
- **Tacky Viscosity** - Rely•a•Bond® Paste is very tacky, which prevents bracket flotation and penetrates the finest mesh bases.
- **Variety of Paste Packaging** - Rely•a•Bond® Paste is available packaged in a jar, preloaded syringe tips or 3.5 g push syringes.
- **Extended Shelf Life (No Refrigeration)** - Rely•a•Bond®, like most Reliance adhesives, has a unique catalyst system that does not require refrigeration.
- **Economical** - Approximately 850 bonds per kit, because only a small amount of primer and paste is required to produce maximum bond strength
- **Suggested Uses** - Bracket adhesive

Rely-a-Bond® Kits

- Rely•a•Bond® kit with Paste in Jar with fluoride..... #RJF
- Rely•a•Bond® kit with Paste in Push Syringes with fluoride..... #RPF
- Rely•a•Bond® kit with Paste in Preloaded Tips with fluoride..... #RTF

Individual Rely-a-Bond® Components

- 14 g Rely•a•Bond® Paste in Jar..... #RPJ
- 14 g Rely•a•Bond® Paste in (4) 3.5 g Push Syringes..... #RPP
- 14 g Rely•a•Bond® Paste in Preloaded Syringe Tips..... #RPT
- 16 g Rely•a•Bond® Primer with fluoride..... #RPRF
- 16 g Rely•a•Bond® Primer without fluoride..... #RPR

Also Available

- 14 g Rely•a•Bond® Thin Paste in (4) 3.5 g Push Syringes (for Indirect Bonding)..... #RPPT

Phase II®

Sealant/Paste System



Contents of Kit

- 18 g Catalyst Paste (A) • 18 g Base Paste (B) •
- 7 g Bonding Resin (A) • 7 g Bonding Resin (B) •
- 9 g Etching Agent (Liquid or Gel) • Mixing Pad •
- Spatulas • Wooden Applicators • Mixing Well •
- Brushes (Sponge pellets may be substituted for brushes) •

Phase II® available with the paste pre-packaged in jars or screw syringes

Product Features

- **Contains Fluoride** - Phase II® Paste and sealant are formulated with a monomer that contains fluoride and will not reduce bond strength.
- **Unsurpassed Bond Strength** - The highly filled paste provides a strong mechanical bond with all types of metal and ceramic bracket bases. All plastic and fiber filled brackets will bond chemically to Phase II® when conditioned with Reliance Plastic Conditioner.
- **Ease of Removal** - Phase II® has a small particle size glass filler which constitutes a high percentage of the material. The paste forms a solid matrix that bonds stronger to itself than it does to the sealant. As a result, the majority of the paste is removed with the bracket and what remains is easier to remove than conventional composites.
- **Smooth, Tacky Consistency** - The pastes mix to a viscosity that is easy to apply to the bracket base and will hold the bracket in position after alignment, without slump. The paste viscosity can be thinned for bonding lingual retainers by adding bonding resin to the mixed paste.
- **No Refrigeration Required** - Phase II®, like most Reliance adhesives, has a unique catalyst system that does not require refrigeration. An active archwire can be placed 5 minutes after last bracket is in place.
- **Suggested Uses** - Bracket adhesive.

Phase II® Kits

Phase II® kit with Paste in Jar with fluoride #PJF

Phase II® kit with Paste in Screw Syringes with fluoride #PSF

NOTE: Assure® or Maximum Cure® Sealants can be substituted for bonding resin in any kit at an additional cost. Please specify.

Individual Phase II® Components

18 g Phase II® Paste A with fluoride in Jar..... #PJAF

18 g Phase II® Paste B in Jar..... #PJB

18 g Phase II® Paste A with fluoride in (3) 6 g Screw Syringes..... #PSAF

18 g Phase II® Paste B in (3) 6 g Screw Syringes..... #PSB

7 g Bonding Resin A with fluoride..... #S7AF

7 g Bonding Resin B..... #S7B

Also Available

20 g Bonding Resin A with fluoride..... #S20AF

20 g Bonding Resin B #S20B

Light Bond™

Light Cure Orthodontic Adhesive



Contents of Kit

- 20 g Paste • (2) 3 cc Sealant Resin •
- 23 g Liquid Etching Agent •
- Mixing Pad • Spatulas •
- Brushes (Sponge pellets may be substituted for brushes) •

Product Features

- **Available With or Without Fluoride** - Light Bond™ Paste and Sealant are formulated with a patented monomer that contains fluoride. It is available without fluoride in Push Syringe only.
- **Unsurpassed Bond Strength** - The Light Bond™ Paste is a highly filled material that delivers maximum bond strength in all situations with all types of appliances
- **Smooth Tacky Paste Viscosity** - The viscosity of Light Bond™ Paste is very tacky which prevents bracket flotation. If a thinner paste viscosity is desired for bonding lingual retainers and acrylic appliances, Light Bond™ Thin or Medium Paste is available. Specify when ordering.
- **Total Paste Polymerization Under All Bracket Bases** - Light Bond™ is comprised of a catalyst and filler that reflects the light to insure complete polymerization under any appliance: metal, ceramic, sapphire or plastic. Once the material has been properly cured, an arch wire can be placed immediately.
- **Ease of Removal** - Light Bond™ is a small particle glass filled material designed specifically for orthodontic use that allows maximum penetration of filler into the finest screen mesh base and removal without enamel scarring.
- **Suggested Uses** - Bracket adhesive

Light Bond™ Kits

Light Bond™ kit with Paste in Push Syringes with fluoride.....#LBEPF

Light Bond™ kit with Paste in Screw Syringes with fluoride.....#LBESF

Light Bond™ kit with Paste in Tips with fluoride #LBETF

NOTE: Light Bond™ Filled Sealant may be substituted for standard Light Bond™ Sealant. Non-fluoride Release Formula available upon request. Light Bond™ Paste is available in Regular, Medium and Thin viscosities. Light Bond™ Paste can be used with any cured sealant. Light Bond™ Kit is available with Assure® Plus All Surface Bonding Resin, Assure® Universal Bonding Resin (hydrophilic) or PRO SEAL® or L.E.D. PRO SEAL® Enamel Protective Sealant as replacements for Light Bond™ Sealants at an additional charge. Blue shade available upon request.

Individual Light Bond™ Components

5 g Light Bond™ Paste in Push Syringe with fluoride#LBPPF

5 g Light Bond™ Paste in Screw Syringe with fluoride.....#LBPSF

5 g Light Bond™ Paste in Tips with fluoride #LBPTF

3 cc Light Bond™ Sealant with fluoride#LBS3F

7 cc Light Bond™ Sealant with fluoride#LBS7F

Pad Lock®

Light Cure Fluorescing Adhesive



Contents of Kit

- 18 g Paste • [2] 3 cc Sealant Resin •
 - 9 g Etching Agent (Liquid or Gel) • Mixing Pad • Spatulas •
 - Brushes (Sponge Pellets may be substituted for Brushes) •
- REFRIGERATION IS RECOMMENDED**

Product Features

- **Bonds Chemically to Metal** - Pad Lock® was developed to increase adhesion to even the finest screen mesh bases. The paste contains a monomer that bonds chemically as well as mechanically to metal bracket bases improving adhesion.
- **Fluoresces** - Pad Lock® Paste will fluoresce when illuminated with any ultraviolet black light. This feature allows the operation to disclose flash not removed at the time of bracket placement and adhesive remnants left at debonding.
- **Single Paste, Light Activated**
- **Compatible With All Sealants and Self Etching Primers**
- **Available With or Without Fluoride**
- **Suggested Uses** - Bracket adhesive.

Pad Lock® Kits

- Pad Lock® kit in Push Syringe with fluoride.....#PLPKF
 Pad Lock® kit in Tips with fluoride.....#PLTKF
 NON-FLUORIDE RELEASE FORMULA AVAILABLE UPON REQUEST.
 NOTE: **Assure® Plus** All Surface Bonding Resin, **Assure®** Universal Bonding Resin, **LED PRO SEAL®** or **PRO SEAL®** Enamel Protective Sealant may be substituted for standard **Light Bond™** Sealant for an additional charge.

Individual Pad Lock® Components

- 4.5 g Pad Lock® Paste in Push Syringe with fluoride.....#PLPF
 4.5 g Pad Lock® Paste in Tips with fluoride.....#PLTF
 3 cc Light Bond™ Sealant with fluoride.....#LBS3F
 3 cc Light Bond™ Filled Sealant with fluoride.....#LBFSF

Also Available

- 7 cc Light Bond™ Sealant with fluoride.....#LBS7F
 7 cc Light Bond™ Filled Sealant with fluoride.....#LBFS7F
 6 cc Assure® Universal Bonding Resin.....#AS
 6 ml PRO SEAL® Light Cure Fluorescing Filled Sealant.....#PRO
 6 ml L.E.D. PRO SEAL® Light Cure Fluorescing Filled Sealant.....#LEDPRO
 UV Disclosing Light.....#UVML

Assure®

Hydrophilic Light Cure Orthodontic Adhesive
Contains Fluoride



Contents of Kit

- 20 g Paste in Screw Syringes •
- [1] 6 cc **Assure®** Universal Bonding Resin •
- 9 g Etchant (Liquid or Gel) • Mixing Pad • Spatulas • Mixing Well •
- Brushes (Sponge Pellets may be substituted for Brushes) •

REFRIGERATION REQUIRED

Product Features

- **Hydrophilic** - **Assure®** Paste is a light cure paste that will bond to slightly contaminated enamel. **Assure®** Paste is the ideal viscosity to prevent bracket flotation and contains fluoride.
- **Bonds to Any Enamel Type** - **Assure®** Universal Bonding Resin will bond to any enamel type, wet or dry, normal or atypical, including deciduous, bleached and dentin. In addition, **Assure®** will bond to composite and metal surfaces without any additional primer.
- **Assure® Universal Bonding Resin** - Can be used under any chemical, light or dual cure adhesive, regardless of manufacturer.
- **Suggested Uses** - Bracket adhesive, Wet field bonding.

Assure® Kit

Assure® kit with fluoride #AKP

Individual Assure® Components

5 g **Assure®** Paste in Screw Syringe with fluoride #APP
6 cc **Assure®** Universal Bonding Resin..... #AS
6 cc **Assure® Plus** All Surface Bonding Resin..... #PLUS

(**Assure®** Universal Bonding Resin may be substituted in kit for **Assure® Plus** All Surface Bonding Resin for an additional cost)

Quick Cure™

Light Cure Bracket Bonding System



Contents of Kit

- 20 g Paste • [2] 3 cc **Light Bond™** Sealant •
- 23 g Liquid Etchant •
- Mixing Pad • Spatulas • Mixing Well •
- Brushes (Sponge Pellets may be substituted for Brushes) •

Product Features

- **Reduced Curing Time** - **Quick Cure™** is a new generation of light cure paste with additional catalyst for a faster and more complete cure with any curing light.
- **Choice of Sealant** - **Quick Cure™** can be used with any light cure sealant and is available in kit form with **Light Bond™** Sealant, **PRO SEAL®** or **Assure®** Universal Bonding Resin for an additional charge.
- **Contains Fluoride** - **Quick Cure™** Paste is formulated with a patented monomer which releases fluoride. University studies have shown that fluoride is continuously released throughout the mouth without reducing bond strength.
- **Suggested Uses** - Bracket adhesive.

Quick Cure™ Kits

Quick Cure™ kit in Push Syringe with **Light Bond™** Sealant..... #QCPLB
Quick Cure™ kit in Push Syringe with **Assure®** Universal Bonding Resin... #QCPAS
Quick Cure™ kit in Tips with **Light Bond™** Sealant..... #QCTLB
Quick Cure™ kit in Tips with **Assure®** Universal Bonding Resin... #QCTAS

NOTE: **Assure® Plus** All Surface Bonding Resin, **Assure®** Universal Bonding Resin, **LED PRO SEAL®** or **PRO SEAL®** Enamel Protective Sealant may be substituted for standard **Light Bond™** Sealant for an additional charge.

Individual Quick Cure™ Components

5 g **Quick Cure™** Paste in Push Syringe with fluoride #QCPP
5 g **Quick Cure™** Paste in Tips with fluoride #QCPT

Protective Eyewear



Uvex Skyper®



Uvex Ultra Spec® 2000

Product Features

More and more orthodontists are using light cured adhesives and newer, more powerful curing lights. It makes sense for you, your assistants and patients to protect your eyes from the powerful light sources now being used everyday in your clinic. Reliance has sourced a very high quality, yet value priced protective eyewear from Uvex. They absorb 99.9% of harmful UVA and UVB radiation, and are available in two designs.

Order Number

Uvex Skyper® - for use without prescription glasses.....#SG

Uvex Ultra Spec® 2000 -
sized to be used over your prescription glasses #SGU

Power Slot™

Quick Curing Fiber Optic Light Guide



Product Features

Now you can significantly increase the power of your curing light by simply replacing the light director. The **Power Slot™** will polymerize light cure bracket adhesives in 10 SECONDS. Reliance **Quick Cure™** Paste can be cured in 6 SECONDS. The **Power Slot™** fits most popular curing lights and has a patented unique flat tip design that maximizes light energy beneath the bracket base. Maximum tapered optics increase the output of your light 2.5 times to insure fast curing of any adhesive with maximum strength.

Order Number

Power Slot™#(Call for corresponding product code)

UV Disclosing Light



Product Features

Take advantage of one of **Pad Lock®** and **PRO SEAL's®** unique features.

PRO SEAL® and **Pad Lock®** will fluoresce under UV light, allowing you to check your application and coverage.

Order Number

UV Disclosing Light #UVML

L.E.D. PRO SEAL® or PRO SEAL®

Light Cure Fluorescing Filled Sealant
Contains Fluoride



Product Features

- **PRO SEAL®** and **L.E.D. PRO SEAL®** are the first sealants that will completely set without an oxygen inhibited layer. This creates a smooth hard surface that protects the enamel.
- **PRO SEAL®** and **L.E.D. PRO SEAL®** are highly filled, no mix light cure sealants that resist toothbrush abrasion.
- **PRO SEAL®** and **L.E.D. PRO SEAL®** also contain a UV fluorescing agent for easy monitoring of sealant coverage.
- **PRO SEAL®** and **L.E.D. PRO SEAL®** can be used under any light cure, chemical cure or dual cure paste system. U.C.S.F. and M.C.G. tests prove the long term benefits of this enamel protective sealant.

Order Number

6 ml **L.E.D. PRO SEAL®**#**LEDPRO**
6 ml **PRO SEAL®**#**PRO**

Maximum Cure® Sealant

Enamel Protective Chemical Cure Sealant System
Contains Fluoride



Product Features

- The traditional self cure bonding agent does not totally polymerize due to oxygen inhibition. Consequently, it is rinsed away at the conclusion of the appliance bonding procedure leaving exposed etched enamel unprotected.
- **Maximum Cure®** Sealant will set in a thin layer. A unique catalyst system prevents oxygen interference during polymerization. The whole labial or lingual tooth surface can be sealed.
- **Maximum Cure®** Sealant is formulated with a patented fluoride monomer.
- Also available: **Maximum Cure®** Filled Sealant for longer lasting protection.
- **CHEMICAL CURE SYSTEM - NO CURING LIGHT REQUIRED:** With **Maximum Cure®** Sealant, multiple teeth can be sealed in the time it would take to cure one tooth with a light cure material.

Order Number

7 g **Maximum Cure®** Sealant A with fluoride#**MC7AF**
7 g **Maximum Cure®** Sealant B#**MC7B**
7 g **Maximum Cure®** Filled Sealant A with fluoride.....#**MCF7AF**
7 g **Maximum Cure®** Filled Sealant B#**MCF7B**
Economy Size **Maximum Cure®** Sealant
20 g **Maximum Cure®** Sealant A with fluoride#**MC20AF**
20 g **Maximum Cure®** Sealant B#**MC20B**

REFRIGERATION REQUIRED

L.C.R.™

Light Cure Retainer Adhesive

**Contents of Kit**

- 20 g Paste • [2] 3 cc Sealant Resin •
- 23 g Liquid Etching Agent • Mixing Pad • Spatulas •
- Brushes [Sponge Pellets may be substituted for Brushes] •

1.5 g L.C.R.™ Paste in Push Syringe with Tips

**Product Features**

- **Longevity and Durability** - L.C.R.™ Light Cured Bonding Paste is a hybrid composite that cures to a smooth glass finish for patient comfort and reduced wear. This highly filled light cured resin paste is an excellent choice when longevity and durability are required. Simply apply according to instructions and light cure.

L.C.R.™ Kits

L.C.R.™ kit with Paste in Push Syringes.....#LCREP

L.C.R.™ kit with Paste in Tips.....#LCRET

NOTE: **Assure® Plus** All Surface Bonding Resin, **Assure®** Universal Bonding Resin, **LED PRO SEAL®** or **PRO SEAL®** Enamel Protective Sealant may be substituted for standard **Light Bond™** Sealant for an additional charge.

Individual L.C.R.™ Components

5 g L.C.R.™ Paste in Push Syringe.....#LCR

5 g L.C.R.™ Paste in Tips.....#LCRT5

Also Available

1.5 g L.C.R.™ Paste in Push Syringe with Tips.....#LCR15

Flow Tain™ and Flow Tain™ L.V.

Light Cure Flowable Microfil Composite

**Suggested Uses**

- Lingual retainer adhesive
- Attachments on thermoplastic tooth aligning products
- Splinting material (e.g. fibers)
- Indirect bonding adhesive
- Arch wire stops

Product Features

- **Flow Tain™** and **Flow Tain™ L.V.** are flowable light cured composites. These materials are very smooth when cured which will reduce wear and enhance patient comfort.
- **Cured Composite Flexibility** - **Flow Tain™** and **Flow Tain™ L.V.**, when cured, have a low modulus of elasticity and are not as brittle as conventional composites.
- **Flowable Viscosity** - Even though **Flow Tain™** and **Flow Tain™ L.V.** are flowable materials, their thixotropic properties prevent them from slumping for better control.
- **Needle Tip Dispenser** - Using the tips provided, these materials are easily injected onto lingual surfaces and across fiber splints.

Individual Flow Tain™ Components1.5 g **Flow Tain™** with Tips.....#FT1.5 g **Flow Tain™ L.V.** with Tips.....#FTLV**Flow Tain™** 19 gauge Tips.....#FTTIPS**Flow Tain™ L.V.** 22 gauge Tips.....#FTLVIPS

Excel®

Specialty Bonding Adhesive
Contains Fluoride



Excel® Regular Kit in Screw Syringes



Excel® Regular Kit in Jars

Contents of Excel® Regular Kit

- 12 g Catalyst Paste (A) • 12 g Base Paste (B) •
- 7 g Bonding Resin (A) • 7 g Bonding Resin (B) •
- 9 g Etching Agent (liquid or gel) • Mixing Pad • Spatulas • Mixing Well •
- Brushes (Sponge pellets may be substituted for brushes) •

Contents of Excel® Fast Kit

- 12 g Catalyst Paste (A) • 12 g Base Paste (B) •
- 9 g Etching Agent (liquid or gel) • Mixing Pad • Spatulas •
- Brushes (Sponge pellets may be substituted for brushes) •

NOTE: Excel® Regular Herbst (for Acrylic Appliances and Splints) contains Reliance Plastic Conditioner for chemical adhesion to all acrylic appliances. Maximum Cure® Sealants can be substituted for sealants in any kit at no additional cost. Excel® Regular Paste is also available in Blue shade. Please specify.

Product Features

- **Contains Fluoride** - Excel® Paste and Bonding Resin are formulated with a patented fluoride releasing monomer that will not reduce bond strength.
- **Unsurpassed Bond Strength** - Excel® is a highly filled adhesive that provides maximum bond strength in a thin film.
- **Thin, Creamy Consistency** - Excel® Regular Set Paste is formulated in a thin viscosity to flow easily into any acrylic appliance, indirect tray or lingual retainer.
- **Extended Working Time** - Excel® Regular Set (3 minute working time) allows the operator sufficient time to spatulate the material, apply it to the appliance, place the appliance into the mouth and clean up excess flash. For individual bracket placement, Excel® Fast Set is recommended, because of it's 60 second working time.
- **Suggested Uses** - Indirect Bonding adhesive, Bonded Acrylic appliances

Excel® Kits

- Excel® Regular kit with Paste in Jar with fluoride.....#ERJF
- Excel® Regular kit with Paste in Screw Syringes with fluoride.....#ERSF
- Excel® Herbst kit for Large Acrylic Appliances and Splints with Paste in Jar with fluoride#EHJF
- Excel® Herbst kit for Large Acrylic Appliances and Splints with Paste in Screw Syringes with fluoride#EHSF
- Excel® Fast kit with Paste in Jar with fluoride#EFJF
- Excel® Fast kit with Paste in Screw Syringes with fluoride#EFSF

Individual Excel® Components

- 7 g Bonding Resin A with fluoride#S7AF
- 7 g Bonding Resin B#S7B
- 12 g Excel® Regular Paste A in Jar with fluoride.....#ERJAF
- 12 g Excel® Regular Paste B in Jar#ERJB
- 12 g Excel® Regular Paste A in (2) 6 g Screw Syringes with fluoride...#ERSAF
- 12 g Excel® Regular Paste B in (2) 6 g Screw Syringes#ERSB
- 12 g Excel® Fast Paste A in Jar with fluoride.....#EFJAF
- 12 g Excel® Fast Paste B in Jar.....#EFJB
- 12 g Excel® Fast Paste A in (2) 6 g Screw Syringes with fluoride....#EFSAF
- 12 g Excel® Fast Paste B in (2) 6 g Screw Syringes.....#EFSB

Phase II® Dual Cure

Sealant / Paste System
Contains Fluoride



Contents of Kit

- 18 g Catalyst Paste [A] with fluoride • 18 g Base Paste [B] •
- [2] 3 cc Sealant Resin • 9 g Etching Agent (Liquid or Gel) •
- Mixing Pad • Spatulas • Mixing Well • Wooden Applicators •
- Brushes [Sponge Pellets may be substituted for Brushes] •

Product Features

- **Combines the Best of Both Chemical and Light Cure Polymerization** - The chemical cure working time is 4 minutes to allow several brackets to be positioned from one mix. An initial set can be immediately achieved by light curing from any angle for 10 seconds. This will prevent the loss of bond strength due to saliva contamination. This unique dual cure process allows complete polymerization under metal brackets where light does not reach.
- **Contains Fluoride** - Phase II® Dual Cure paste and sealant are formulated with a patented fluoride releasing monomer. Phase II® Dual Cure can be used with any chemical or light cure sealant.
- **Suggested Uses** - Bracket adhesive.

Phase II® Dual Cure Kits

Phase II® Dual Cure kit with Paste in Jar with fluoride #DCJF

Phase II® Dual Cure kit with Paste in Screw Syringes with fluoride.... #DCSF

NOTE: **Assure® Plus** All Surface Bonding Resin, **Assure®** Universal Bonding Resin, **LED PRO SEAL®** or **PRO SEAL®** Enamel Protective Sealant may be substituted for standard **Light Bond™** Sealant for an additional charge.

Individual Phase II® Components

18 g Phase II® Dual Cure Paste A in Jar with fluoride..... #DCJAF

18 g Phase II® Paste B in Jar..... #PJB

18 g Phase II® Dual Cure Paste A in [3] 6 g Screw Syringes with fluoride #DCSAF

18 g Phase II® Paste B in [3] 6 g Screw Syringes..... #PSB

3 cc **Light Bond™** Sealant with fluoride #LBS3F

7 g Bonding Resin A with fluoride..... #S7AF

7 g Bonding Resin B..... #S7B

Band-Lok®

Two Paste Compomer Band Cement
Contains Fluoride



Contents of Kit

- 24 g Paste A • 24 g Paste B •
- Mixing Pad • Spatulas •

Product Features

- **No Messy Powder-Liquids** - Consistent mix viscosity every time.
- **Bonds to Metal** - Band-Lok® will bond chemically to stainless steel bands and crowns. Removal is all in one piece with cement in the appliance.
- **Extended Working Time** - Four minute working time at room temperature without frozen slab.
- **Dual Curable** - Band-Lok® is a dual cure cement. It can be used with or without dental curing light. After flash is cleaned a 20-second light activation allows full protection from moisture. Without light activation a complete set is achieved in 10 minutes.

Available in Natural or Blue Shade

Band-Lok® Kits

Band-Lok® kit in Natural shade with fluoride#BL
Band-Lok® kit in High Visibility Blue shade with fluoride.....#BLB

Individual Band-Lok® Components

6 g Band-Lok® Paste A in Screw Syringe with fluoride#BLSA
6 g Band-Lok® Paste B in Screw Syringe#BLSB
6 g Band-Lok® Paste A Blue shade in Screw Syringe with fluoride ..#BLSAB

Crown Lok™

Luting Cement for Bands and Metal Crowns



Contents of Kit

- 24 g Catalyst Paste A • 24 g Base Paste B •
 - 9 g Etching Agent (Liquid or Gel) • Mixing Pad • Spatulas •
- REFRIGERATION RECOMMENDED**

Product Features

- **Self Curing** - Crown Lok™ is a self curing, two minute working time fluoride releasing cement. This versatile cement will bond chemically to stainless steel crowns and bands, reducing the chance of failure and facilitation cleanup at the conclusion of treatment

Available in Natural or Blue Shade

Crown Lok™ Kits

Crown Lok™ kit in Natural shade with fluoride#CL
Crown Lok™ kit in High Visibility Blue shade with fluoride.....#CLB

Individual Crown Lok™ Components

6 g Crown Lok™ Paste A in Syringe with fluoride#CLSA
6 g Crown Lok™ Paste B in Syringe#CLSB
6 g Crown Lok™ Paste A Blue Shade in Screw Syringe with fluoride..#CLSAB

Self Mixing Dual Syringes for Band-Lok® and Crown Lok™

Available with Band-Lok®, Band-Lok® Blue, and Crown Lok™ Blue paste



Contents of Band-Lok® Kit

- (4) Dual Syringes of Paste • (60) Mixing Tips •

Contents of Crown Lok™ Blue Kit

- (4) Dual Syringes of Paste • (60) Mixing Tips •
- 9 g Liquid Etch •

Product Features


- These syringes mix and dispense the paste in one easy step! The A & B syringes are joined and use special mixing canula tips, which eliminate operator mixing and provide a consistent, automatic mix every time!

Self Mixing Dual Syringe Kits

Band-Lok® Natural shade with fluoride#SMBLK
Band-Lok® Blue shade with fluoride#SMBLBK
Crown Lok™ Blue shade with fluoride#SMCLBK

Ultra Band-Lok®

Single Paste Compomer Band Cement
Contains Fluoride

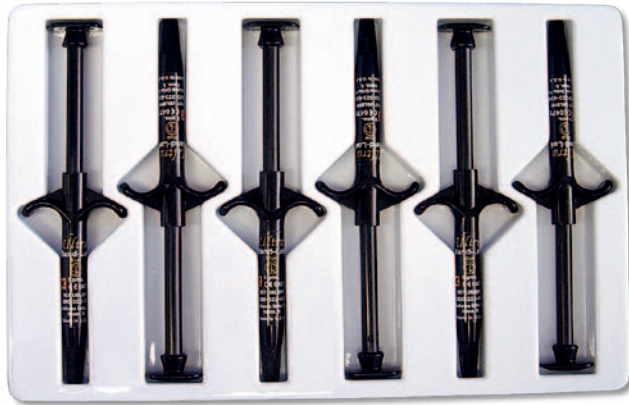

Also available in Luer Loc Syringe



Contents of Ultra Band-Lok® Kit

- 30 g Paste in Push Syringes •

Contents of Luer Loc Ultra Band-Lok® Kit

- 30 g Paste in Luer Loc Push Syringes •
- (60) 15 gauge Luer Loc Dispensing Tips •

Product Features

- **Single Paste System** - No wasteful mixing
- **Light Activated** - No working time constraints
- **No Mixing**
- **All Purpose Cement** - Ultra Band-Lok® can also be used to bond large acrylic appliances and as an occlusal build-up material. When bonding the appliances you must first etch the enamel surface.
- **Bonds to Metal** - Ultra Band-Lok® will bond chemically to stainless steel bands. Removal is all in one piece with cement in the appliance

Available in Natural or Blue Shade

Ultra Band-Lok® Kits

- Ultra Band-Lok® kit in Natural shade with fluoride.....#UBL
- Ultra Band-Lok® kit in Hi Vis Blue shade with fluoride.....#UBLB
- Ultra Band-Lok® Luer Loc w/Dispensing Tips with fluoride.....#LLUBL
- Ultra Band-Lok® Luer Loc Hi Vis Blue shade w/Dispensing Tips with fluoride.....#LLUBLB

Individual Ultra Band-Lok® Components

- 5 g Ultra Band-Lok® Syringe with fluoride.....#UBLPP
- 5 g Ultra Band-Lok® Hi Vis Blue shade Syringe with fluoride..#UBLBPP
- 5 g Ultra Band-Lok® Luer Loc Syringe w/Dispensing Tips with fluoride..#LLUBLP
- 5 g Ultra Band-Lok® Luer Loc Hi Vis Blue shade Syringe w/Dispensing Tips with fluoride.....#LLUBLBP
- (60) 15 gauge Luer Loc Dispensing Tips.....#LLTIPS



Starter Kit



Economy Kit

Precedent®

Glass Ionomer Cement
Contains Fluoride

Product Features

- **High Compressive & Tensile Strength** - Unlike zinc phosphate cements, Precedent® will not wash out from under orthodontic bands
- **Reduce Likelihood of Decalcification** - Precedent® cement leaches fluoride ions to reduce decalcification.
- **Bonds in a Wet or Dry Environment** - Precedent® is a hydrophilic glass ionomer that will bond in a saliva or water contaminated field.
- **Suggested Uses** - Band cement

Precedent® Kits

- Precedent® Starter kit.....#PRS
- Precedent® Economy kit.....#PRE

Individual Precedent® Components

- 15 g Precedent® Powder.....#PRP15
- 100 g Precedent® Powder.....#PRP100
- 15 cc Precedent® Liquid.....#PRL15
- 90 cc Precedent® Liquid.....#PRL90

Contents of Kit

Starter Kit

- 15 g Powder •
- 15 cc Liquid •
- Measuring Scoop •
- Spatulas and Mixing Pad •

Economy Kit

- 100 g Powder •
- 90 cc Liquid •
- Measuring Scoop •
- Spatulas and Mixing Pad •

O.V.S.™

Orthodontic Veneer System



Before OVS Treatment



After OVS Treatment

Contents of Kit

- **Faci Form** Crown Form #12 Left (5 med & 5 large) •
- **Faci Form** Crown Form #22 Right (5 med & 5 large) •
- 2.5 g **Ortho Veneer™** Paste Shade A-1 • 2.5 g **Ortho Veneer™** Paste Shade A-2 •
- 2.5 g **Ortho Veneer™** Paste Shade B-1 •
- 3 cc **Light Bond™** Sealant • (16) **Stick-n-Place** instruments •

Product Features

- **An Excellent Solution to a Unique Problem!** - O.V.S.™ has been perfected and used by Dr. James Hilgers and Dr. Debra Cook in their own practice. A common problem when treating orthodontic cases is the need to bond an attachment to a peg lateral or an undersized tooth. With O.V.S.™ that problem has been solved! O.V.S.™ will provide a complete in-office system that creates an aesthetic and properly sized temporary veneer surface to bond your attachment. This saves you and your patients both time and money by eliminating the need for a full temporary crown.

O.V.S.™ Kit and Individual Components

O.V.S.™ kit.....	#OVS
Faci Form™ Crown Form #12 Left (5 med & 5 large).....	#FFL
Faci Form™ Crown Form #22 Right (5 med & 5 large).....	#FFR
2.5 g Ortho Veneer™ Paste Shade A-1	#OVPA1
2.5 g Ortho Veneer™ Paste Shade A-2	#OVPA2
2.5 g Ortho Veneer™ Paste Shade B-1	#OVPB1
3 cc Light Bond™ Sealant with fluoride.....	#LBS3F
(16) Stick-n-Place™ instruments	#SNP20

Therma-Cure™

Heat Curable Lab Adhesive for Custom Base Indirect Bonding



Product Features

- **Tacky Paste Viscosity** - No bracket flotation.
- **Unlimited Working Time** - For bracket positioning and excess flash cleanup.
- **Heat Curable** - Simply place in toaster oven for 15 minutes at 325°F.
- **Highly Filled Composite** - Produces maximum bond strength with all metal, ceramic and plastic brackets. Simply bond in mouth with sealant.

Order Number

5 g Therma-Cure™ syringe #TC

Custom I.Q.™

Fluoride Releasing Fast Set Indirect Bonding Sealant



Product Features

- **Easy Application** - Custom I.Q.™ is an A & B sealant that is applied to the prepared tooth and custom bracket base without mixing. Simply apply Part A to the custom bracket base and part B to the enamel and seat tray.
- **Fast Set Time** - The indirect tray is held in position for one minute and then passively rests for four additional minutes.
- **Chemical Cure System** - No curing lights required! Bonds to any custom base pad when used with Reliance **Plastic Conditioner** or **Enhance™ Adhesion Booster**.

Refrigeration Required

Order Number

7 g Custom I.Q.™ Sealant A with fluoride #CIQA

7 g Custom I.Q.™ Sealant B #CIQB

Adhesion Boosters

Enhance™

Adhesion Booster for Atypical or Contaminated Enamel



Product Features

- **Enhance™ Adhesion Booster** will increase the adhesion of any chemical, light or dual cure system to the following surfaces: fluorosed, hypo-calcified or deciduous enamel, acrylic, bleached or wet enamel, dentin, metal surfaces, composite restorations and enamel.

Order Number

3 cc Enhance™ Adhesion Booster Part A #EBAA

3 cc Enhance™ Adhesion Booster Part B #EBAB

6 cc Enhance™ Adhesion Booster Part A #EBAA6

6 cc Enhance™ Adhesion Booster Part B #EBAB6

Enhance™ L.C.

Single Component Enhance Light Cure



Product Features

- Provides all the benefits of chemical cure **Enhance™** without the mixing. Must be used with light cure sealant. This single liquid (no mixing) is applied to the etched enamel surface then lightly air dried. **Light Bond™ Sealant** is placed directly onto the **Enhance™ L.C.** and light cured for 10 seconds (20 seconds required if **PRO SEAL®** is used).
- Most convenient for **Light Bond™** and **Pad Lock®** users!

Order Number

6 cc Enhance™ Light Cure #LCE6

Assure® Plus

All Surface Bonding Resin



We are proud to release our newest version of **Assure®** that will simplify the primer market – **Assure® Plus**. In addition to all the substrates our original **Assure®** bonds to, **Assure® Plus** will bond to zirconia and acrylic temps without any additional primers and to porcelain without hydrofluoric acid.

Product Features

• ONE BOTTLE FOR ALL SURFACES

- Wet or Dry Normal Enamel
- Wet or Dry Atypical Enamel
- Dentin / Cementum
- Composite
- Gold, Amalgam, Stainless Steel
- Zirconia
- Acrylic Temporary / Pontic Teeth
- Porcelain

Order Number

6 cc **Assure® Plus** All Surface Bonding Resin.....#PLUS

2 ml **Assure® Plus** All Surface Bonding Resin#PLUS2ML

Assure®

Universal Bonding Resin



Product Features

- **Assure®** Universal Bonding Resin will increase adhesion to normal, atypical, dry or slightly contaminated enamel. This universal bonding agent bonds to metal restorations without additional primers. No drying or light curing required and can be used with any light or chemical cure system.
- **Assure®** Universal Bonding Resin will bond to the following surfaces without additional primers: fluorosed, deciduous, bleached or wet enamel, dentin, metal surfaces and composite restorations.
- **Assure®** Universal Bonding Resin is **Enhance™** Adhesion Booster and bonding resin in one step!

Order Number

6 cc **Assure®** Universal Bonding Resin.....#AS

2 ml **Assure®** Universal Bonding Resin..... #AS2ML

Reliance S.E.P.™

Self Etching Primer



4 ml (mixed etchant and primer) #SEPPRO

Blue Gel Etchant In Syringe

High visibility, high viscosity blue gel enamel etchant, with anti-microbial to prevent bacteria buildup after etching



Contains 25 g Gel Etchant in (5) Syringes with 50 tips #BEGE

Etchant Liquid

9 g Etchant Liquid #EL9

Economy Etchant Liquid23 g Economy Etchant Liquid #EL24
23 g Cherry Flavor Economy Etchant Liquid #EL24C**Etchant Gel**

9 g Etchant Gel #EG9

Economy Etchant Gel

18 g Economy Etchant Gel #EG18

Gel Etchant In Syringe

Contains 32.5 g Gel Etchant in (5) Syringes with 100 tips #EGES

Surface Conditioners

Porc-Etch™

Porcelain Etching Agent with Barrier Gel



7 g Porc-Etch™ #POE

Porcelain Conditioner

For bonding to Ceramic Crowns



15 cc Porcelain Conditioner... #PO

Rescue™

All Surface Bonding Kit

Rescue is the answer for bonding to enamel, composite, porcelain, plastic and metal surfaces... all in one kit!

**Contents of Rescue All Surface Bonding Kit**

- 15 cc Porcelain Conditioner •
- 7 g Porc Etch with Barrier Gel •
- 10 cc Plastic Conditioner •
- 6 cc Assure® Universal Bonding Resin •

Rescue All Surface Bonding Kit..... #R

Plastic Conditioner

For bonding to Composite Restorations and Acrylic Appliances



10 cc Plastic Conditioner..... #PL

Metal Primer

For bonding to Gold, Stainless Steel and Amalgam



8 ml Metal Primer #ME

1st & Final®

Pre-Packaged Orthodontic Pumice



Product Features

- **Pre-Packaged** - 200 sterile cups for individual patients use. Pre-mixed – no messy mixing with water. Just open and apply with prophyl cup.
- **Improved Bond Strength** - University test shows that 1st & Final® will pumice more effectively and rinse cleaner for increased bond strength.
- **No Added Chemicals** - No flavoring, fluorides or oils which could leave a film that may inhibit or weaken the bond.
- **Enamel Polisher** - 1st & Final® can be used as the final polishing step with the Reliance Renew™ Finishing System.

Order Number

1st & Final® 200 Sterile Cups.....#FF

Prep Step™

Disposable Propyl Angle



Product Features

- 144 individually wrapped and disposable Prep Step™ Propyl Angles are manufactured to run smooth and vibration free. Latex free, 4 web and soft cup are designed to fit a slow speed hand piece. Just throw away after use!

Order Number

Prep Step™ 144 Propyl Angles.....#DPA



Access™

Cheek Retractor



Product Features

- **Unparalleled Access** - Special posterior pads comfortably hold the cheeks away and allow excellent accessibility to the bicuspid and molar area.
- **Tension adjustment** - Spring stainless steel bows allow for adjustment and pivot points provide greater comfort.
- **New Generation Plastic** - The Access™ cheek retractor uses a proprietary plastic and stainless steel bow to assure unsurpassed durability and dependable use over and over again. Separates into three parts for easy cleaning. The Access™ is a simple and elegantly engineered plastic cheek retractor that can be chemically or heat sterilized in an autoclave or dry heat sterilized! Adult and Pedodontic sizes assure the proper fit for any patient.

Order Number

Access™ Adult Cheek Retractor.....#CH1
Access™ Pedo Cheek Retractor.....#CH2

Renew® Finishing System Burs

These Burs are the **FIRST STEP** in the **Renew® Finishing System**

Renew® Bur 118S



- The **118S** is a special 18 fluted bur design for controlled removal of bulk bonding adhesive after the brackets have been removed. These non-tapered friction grip burs cut quickly and do not promote fissures.

Renew® Bur 118S.....#RSB

Renew® Bur 218



- The **218** bur is an 18 fluted design that is tapered to reach into those tight areas along the gingival and occlusal when removing residual adhesive. The wide flutes self clean.

Renew® Bur 218.....#RSB2

Renew® Bur 118L



- The **118L** is a longer shank version of the **118S** with the same excellent adhesive removal properties. These non-tapered friction grip burs cut quickly and do not promote fissures.

Renew® Bur 118L.....#RSBL

Renew® Lingual Bur 815



- The **815** Lingual composite bur with its friction grip, custom length and cutting head, make this the ideal carbide bur for lingual composite removal.

Renew® Lingual Bur 815.....#RSB815

Renew® Finishing System Points #383

The #383 Point is the **SECOND STEP** in the **Renew® Finishing System**

- The **Renew® Finishing System** Points will remove small residual amounts from the tooth surface and will polish the enamel without damage. It's rubberized, tapered shape facilitates removal of adhesive close to the gingival. The #383 Point is the **SECOND STEP** in the **Renew® Finishing System**. The #383 points are available in Straight Hand piece, Latch or Friction Grip. Each pack contains 12 points. Please specify.

Straight Hand Piece



Straight Hand Piece pack of 12#RSPHP

Latch



Latch pack of 12#RSPL

Friction Grip



Friction Grip pack of 12#RSPFG

Diamond Bur 129

- A medium abrasion diamond bur recommended for roughening and preparing plastic, composite restorations and metal surfaces. Reliance now offers an excellent quality medium diamond bur at a down to earth price.



Diamond Bur 129.....#DB


Restore™

Diamond Polishing Paste for metal, composite and porcelain restoration



2 g Syringe Restore™#RESPP2

The Complete Orthodontic Bonding Reference Guide

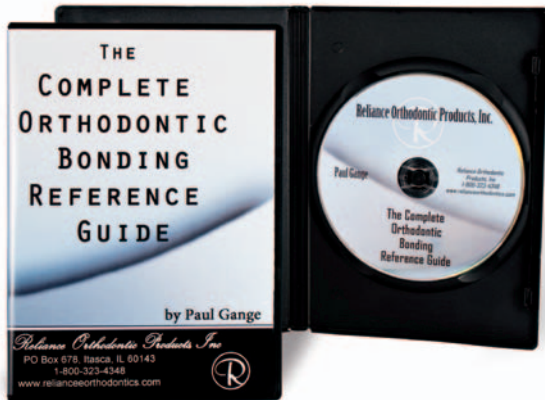

Product Features

- Developed by Paul Gange, respected expert on orthodontic bonding. The Complete Orthodontic Bonding Reference Guide covers every bonding topic and includes tooth preparation, bonding, debonding, lingual retention, indirect bonding and everything in between! Ideal for training chair side staff. Excellent reference for doctor and experienced assistants. Comprehensive, succinct and easy to use, this state-of-the-art reference guide is available in DVD format for your convenience.

Reference Guide Index

- Section 1 – Enamel Preparation
- Section 2 – Direct Bonding Adhesives
- Section 3 – Bonding Troubleshooting
- Section 4 – Bonding to Composite Restorations
- Section 5 – Bonding to Metal Surfaces
- Section 6 – Bonding to Porcelain
- Section 7 – Rebonding using the Same Bracket
- Section 8 – Treatment of Undersized Teeth:
The Orthodontic Veneer System
- Section 9 – Bonding in a Wet Field
- Section 10 – Bonding Lingual Retainers and Bite Turbos
- Section 11 – Band Cementation
- Section 12 – Crown Cementation
- Section 13 – Bonding Large Acrylic Appliances
- Section 14 – Aligner Attachments
- Section 15 – Debonding and Adhesive Removal
- Section 16 – Indirect Bonding

Order Number

Orthodontic Bonding Reference Guide DVD Format..... #DOBRG2

THIS ITEM IS NON-RETURNABLE.

Chairside Bonding Technique Cards



Product Features

- Informative and convenient – all bonding procedures at your fingertips! Every possible intra-oral bonding situation is covered. These durable laminated cards are attached on a ring so they can be kept chair side.

Order Number

Chairside Bonding Technique Cards #TECH

THIS ITEM IS NON-RETURNABLE.



RETAINIUM®

The Next Generation Flat Metal Lingual Retainer
Patent Pending

**Product Features**

- **Single Strand** - 70% stronger than multi strand braided wires.
- **Flat Construction** - maintains torque control without interference.
- **Nickel Free Titanium** - Safe.
- **Uses:** splinting of upper or lower incisors, hold diastemas, maintain extraction sites, periodontal splinting of compromised teeth

Order Number

RETAINIUM™ 10 (6") lingual retainer wires #RET

Stainless Steel Ortho FlexTech® Dental Retention System

The Gold Standard without the Gold Price

**Product Features**

- 30% increase in tensile strength.*
- 300% increase in torsional strength.*
- Half the cost of GOLD Ortho-Flextech®.
- Same great patented design and ease of use for superior results.

* Compared to gold Ortho-FlexTech®.

Order Number

Stainless Steel Ortho FlexTech® - 30" stainless steel chain . #SSOFT

Bond•a•Braid®

Lingual Retainer Wire

**Order Number**

Bond•a•Braid® 10 (6") lingual retainer wires #BAB2

Ortho FlexTech® Dental Retention System

The Next Generation of Lingual Retention Using Solid White Gold Chain

**Order Number**

Ortho FlexTech® - 30" gold chain..... #OFT2

Kobayashi Tie Hooks**Order Number**

Kobayashi Tie Hooks - .010 100/Pkg. #K10
Kobayashi Tie Hooks - .012 100/Pkg. #K12
Kobayashi Tie Hooks - .014 100/Pkg. #K14

Pre-Twisted Kobayashi Tie Hooks

Wires are chemically treated and spot welded together
to prevent breakage

**Order Number**

Pre-Twisted Kobayashi Tie Hooks - .012 100/Pkg. #SK12
Pre-Twisted Kobayashi Tie Hooks - .014 100/Pkg. #SK14

Disposable Benda Brushes



1 1/2" Benda Brushes (250).....#BB250

Brush Handle And Tips



Brush Handle.....#BH
Brush Tips (250).....#BT250

Sponge Pellets



Sponge Pellets (500).....#SP

Micro Brushes



Micro Brushes (100).....#MB
Super Fine Micro Brushes (100).....#SFMB

Nylon Mixing Spatulas



Spatulas (50)#S

Mixing Pads




Non-Absorbant Mixing Pads (3).....#PADS

Mixing Wells



Mixing Wells (5).....#MXW

S.E.P. Mixing Wells

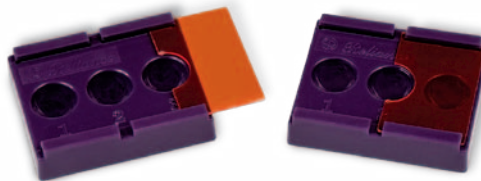

S.E.P. Mixing Well with Light
Impervious Shield (3).....#SEPW

Mark II CR Syringe



Mark II CR Syringe#CRII

Portable Chairside Sandblast Units

Available in a 4-hole adapter which screws directly to your chair side handpiece connector...



- New round shape to improve surface finish.
- Now fabricated entirely out of Stainless Steel thereby increasing its resistance to discoloration and scratching.
- The fit of the gasket around the adapter pins has been changed so silicone is no longer required to secure the gasket to the adapter.

Contents of Kit

- 4 Hole Handpiece Attachment • (5) Regular and (5) Lingual Small BE/50 Micron Hard Abrasive Tips •

4 Hole Handpiece Attachment.....#EMHA

Hard Abrasive Tips



25 Large BE/50µ Hard Abrasive Tips
Regular#EMTRL
Lingual#EMTLL



25 Small BE/50µ Hard Abrasive Tips
Regular#EMTRS
Lingual#EMTSL

Bond Aligner®

Bonds any Metal, plastic or ceramic, regardless of manufacturer
(Patented)



Product Features

A clear light cure paste that will bond any metal, plastic or ceramic button hook or power arm to any aligner material, regardless of manufacturer. In addition, **Bond Aligner** can be used to fabricate composite buttons in a mini mold form.

Order Number

1.4 g **Bond Aligner™** syringe with 5 tips.....#BA

Aligner Pontic Paint

Perfect•A•Smile®

Mask aesthetically compromised teeth during treatment
(Patented)



**Additional Shades
Available Soon!**



Product Features

Single light cure paste will bond to any aligner material to mask diastemas, rotations or aesthetically compromised teeth during treatment. Cured material will not crack or splinter regardless of distortion of the aligner. Available in 3 shades: Light, Medium and Dark.

Contents of Kit

- [1] 8 g Pontic Paint Shade Light •
- [1] 8 g Pontic Paint Shade Medium •
- [1] 8 g Pontic Paint Shade Dark •
- 1 Pack of 50 Brushes • 3 Spatulas •

Perfect•A•Smile™ Kit and Individual Components

Perfect•A•Smile™ Kit	#PAS
Perfect•A•Smile™ shade Light 8 g	#PPB1
Perfect•A•Smile™ shade Medium 8 g	#PPA2
Perfect•A•Smile™ shade Dark 8 g	#PPA3


A²™ Aligner Attachment Kit

Everything you need to fabricate oval or rectangular aligner attachments



Contents of Kit

- 5 Oval Attachment Forms •
- 5 Rectangular Attachment Forms •
- 2 ml **Assure**® Universal Bonding Resin •
- 1.5 g **LCR**™ Paste • Attachment Holder •

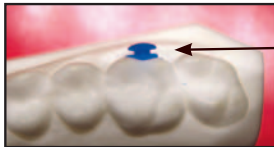

Button Hook Attachment
For Class II Elastics
(sold separately)

Product Features

- **Oval Attachment Form** - Can be used for the retention of thermoplastic aligners
- **Rectangular Attachment Form** - Can be used for individual tooth movement and for the extrusion of teeth
- **Assure**® **Universal Bonding Resin** - Can be used under any chemical, light or dual cure adhesive, regardless of manufacturer, and will bond to all intraoral surfaces
- **LCR**™ - Cures smooth as glass for patient comfort.

A² Aligner Attachment™ Kit

A ² ™ Aligner Attachment kit.....	#A2
5 Oval Attachment Forms.....	#AO
5 Rectangular Attachment Forms.....	#AR
2 ml Assure ® Universal Bonding Resin	#AS2ML
Attachment Holder.....	#BH
1.5 g LCR ™ Paste.....	#LCR15

Also Available

Button Hook Attachment	#ABUT
A ² ™ Disposable Tips (5)	#A2TIP

Oral Surgery Extrusion System

C&C™ Disposable Composite and Chain Kit



Everything you need for a more successful and less painful extrusion system



Contents of Kit

- 2 Disposable Pad Attachment Forms with Chain •
- 3 Disposable Oval Attachment Forms with Chain •
- 2 ml **Assure**® Universal Bonding Resin •
- 1.5 g **FlowTain**™ Flowable Composite •

Product Features

- **More Comfort** - The **C&C**™ Kit does not use protruding metal attachments. Our patent pending design results in more comfort for the patient during treatment AND at the time of removal.
- **Easier Removal** - Instead of "pulling" the attachment off, as with conventional bonded pads, it is removed using a bur, which results in less trauma and discomfort for the patient.
- **Cures Smooth** - Unlike a bracket, the composite cures as smooth as glass and has a very low profile for patient comfort.

C&C™ Kit and Individual Components

C&C™ Disposable Composite and Chain Kit.....	#DCCK
1 Disposable Pad Attachment with Chain.....	#DCCP
1 Disposable Oval Attachment with Chain	#DCCO
2 ml Assure ® Universal Bonding Resin	#AS2ML
1.5 g FlowTain ™ Flowable C&C Paste with Tips	#FTCC

Also Available

C&C™ Disposable Chemical Cure Self Mixing Syringe Kit	#DCCCK
Chemical Cure Self Mixing Syringe and Tips (6g Paste A, 6g Paste B) ..	#SMEF






Reliance Orthodontic Products, Inc.

www.relianceorthodontics.com

(800) 323-4348 • (630) 773-4009 • Fax (630) 250-7704

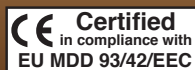

US FDA
REGISTERED

Certificates available upon request

Health Canada
REGISTERED

